

MedCare *ECRAN FACIAL RELEVABLE* RETRACTABLE FACE SHIELD

www.vmedcare.co



Food Services

SAFEGUARD YOUR SMILE



Commercial



Retail

Unlimited applications:



Transportation



Schools



Manufacturing

Retractable Visor

- Visor can be lifted up and stopped at 90 degrees
- Convenient and comfortable for users

Replaceable Lens

- Quick and easy lens replacement
- More eco-friendly
- Kit includes visor and replacement lenses



Made in Canada

Elastic Strap
Long lasting and comfortable fit



Headband

For increased comfort and extended wear


Visor

Design allows space for glasses, goggles and surgical masks (including N95)

Lens

- Anti-fog and anti-static
- Optically clear and distortion-free vision
- Easy replacement
- Lens can pivot 90 degrees

Lift the visor and stop at 90 degrees

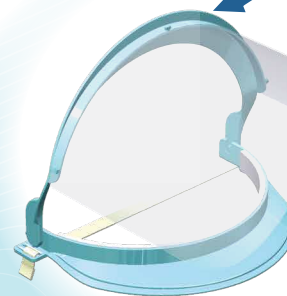


Retractable Visor

3 STEPS OF LENS REPLACEMENT


STEP 1

Remove lens from visor

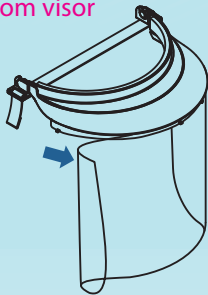


STEP 2

With protective film facing outwards, assemble new lens by securing onto visor

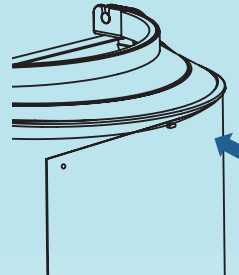

STEP 3

Peel tab to remove protective film

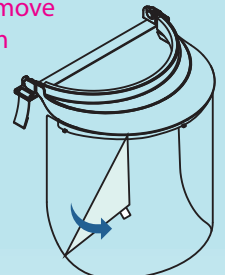

SPECIFICATIONS

Dimensions	Outer Length x Width (mm)	222 x 160
	Inner Length x Width (mm)	160 x 95
	Height (mm)	210
	Lens Width (mm)	345
Weight	Overall (gram)	71
Materials	Headpiece	Polypropylene
	Head band	Polyurethane
	Lens	PET
	Stretch strap	Polyester
Lens Treatment	Anti-fog, Anti-static	